

Push-To-Set™ Pediatric Intermittent Suction Unit, Analog (ISO/CE)

DESCRIPTION AND APPLICATION

The Ohio Medical® Pediatric Push-To-Set™ Intermittent Suction Unit (PED-ISU) is a medical suction device designed to deliver continuous or intermittent controlled vacuum from a hospital pipeline vacuum system. It provides the versatility of an intermittent mode for nasogastric drainage and a continuous mode for tracheal, pharyngeal, or surgical suctioning. The Pediatric model is designed to prevent vacuum adjustment above a preset vacuum limit level. This device is ideal for neonatal and pediatric patients or sensitive applications where limiting the peak vacuum level is important. The intermittent mode uses a pneumatic logic device to turn the suction "ON" and "OFF" automatically (intermittent starts in the "ON" cycle). This unique Unilogic timing module sets the industry standard for quiet, reliable, controlled intermittent suction. The Pediatric model is adjustable from 0-18.0 kPa (0-135 mmHg), to meet most low clinical suction requirements. It also contains a positive pressure relief to reduce potential injury to the patient and equipment if momentarily attached to a positive pressure gas source.

The Pediatric Push-To-Set™ Intermittent Suction Unit is a low vacuum/high flow regulator.

FEATURES AND BENEFITS

- **EXCLUSIVE!** Patented Push-To-Set™
Improves patient safety
May prevent inadvertent unregulated suction
Simple, one handed operation
- **Baby Icon indicates Pediatric or Neonatal Limited Model**
- **Negative Pressure Relief Valve**
Pediatric model is factory preset to 18.0 kPa (135 mmHg)
- **Unilogic Timing Module**
Pneumatic, quiet, reliable operation
Returns to atmospheric pressure during "OFF" cycle
Timing cycles
Independently adjustable 3 to 30 seconds "ON" and "OFF" cycles
Factory preset at 15 seconds on/8 seconds off (± 3 seconds)
Begins in the "ON" mode
- **Flow Rates**
(Intermittent) Preset at 8 L/min (ASTM)
(Continuous) 0 to 80 L/min (without fittings at full increase)
- **Large, Easy to Read Gauge**
Color coded ranges: Low(green), Medium(yellow), High(orange)
Analog: Familiar, industry-standard analog design
- **Safety Features**
Push-to-Set™ automatic occlusion design
Positive pressure safety relief valve
- **Modular Components**
Common components reduce service part inventory
Simplifies service procedures
Low cost replacement parts and kits
- **Durable Design and Construction**
Few mechanical parts
Rugged, shatter-resistant case
Corrosion & Lubrication Free
Service-free backplate
- **CE marked**
- **Ordering Options**
Available in various configurations to including adapters for connection to source vacuum. Also available with options for collection devices and overflow safety protection.



SPECIFICATIONS*

Flow Rate

Intermittent – Preset at 8 L/min per ASTM
 Continuous – 0-80 L/min without fittings at full increase

Timing

3 to 30 seconds
 Independent on and off
 Factory settings, Intermittent timing: 15 seconds on/8 seconds off (± 3 seconds)
 Begins in the “ON” mode

Vacuum Range

0-18.0 kPa (0-135 mmHg)

Gauge Range

0-21.3 kPa (0-160 mmHg)

Gauge Accuracy

±5% Full Scale Deflection or ± 1.3 kPa (8 mmHg)

Dimensions

16.5 x 7.1 x 12.2 cm
 (6.5” H X 2.8” W X 4.8” D)

Weight

0.57 Kg (1 lb 4oz) less fittings

Applied Standards

ASTM F960 and ISO 10079-3

* Specifications are nominal, subject to change without notice. Please see “Instructions for Use” manual for the most current regulated vacuum ranges.

SHIPPING INFORMATION

Shipping Weight

0.77 Kg (1 lb 11 oz) less fittings and adapters

Package Size


16.5 x 26.67 x 20.32 cm
 (6.5” H x 10.5” W x 8” D)

Packaging

One per carton



Parts Configuration

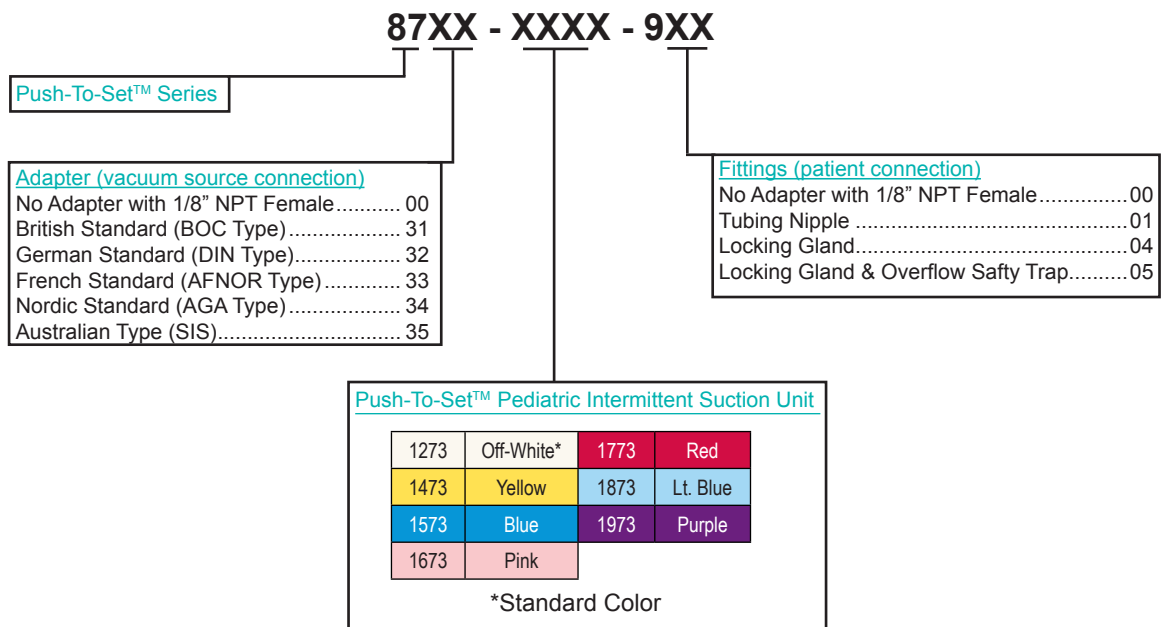

© 2017 Ohio Medical LLC. All rights reserved.
 This document contains information that is proprietary and confidential to Ohio Medical LLC.
 Use of this information is under license from Ohio Medical LLC. Any use other than that authorized by Ohio Medical LLC is prohibited.



Ohio Medical and the Ohio Medical logo are registered trademarks and Push-To-Set is a trademark of Ohio Medical LLC. U.S. Patent #6,623,463 (PushTo-Set™ Vacuum System).

Ohio Medical Authorized Representative

