

Intermittent Suction Unit Adult

DESCRIPTION AND APPLICATION

The Ohio Medical® Intermittent Suction Unit (ISU) is a combination continuous and intermittent regulator for use with hospital pipeline vacuum systems. It provides the versatility of an intermittent mode for nasogastric drainage and a continuous mode for tracheal, pharyngeal, or surgical suctioning. In the intermittent mode, the unit uses a pneumatic logic device to alternately turn the suction “ON” and “OFF.” This unique Unilogic Module provides quiet and reliable controlled intermittent suction. The unit is equipped with independent controls to adjust the duration of the “ON” and “OFF” cycles, the level of vacuum, and the rate of flow.

The regulator also provides continuous vacuum, adjustable over the full range of available vacuum, to meet all clinical continuous suction requirements.

The Intermittent Suction Unit is classified as a high vacuum/high flow regulator.

FEATURES AND BENEFITS

- **Unilogic Module**
 - Pneumatic
 - Quiet, reliable operation
 - Returns to atmospheric pressure during “OFF” cycle
- **Fully Adjustable**
 - Timing cycles
 - 3 to 30 seconds
 - “ON” and “OFF” independently adjustable
 - Factory preset at 15 seconds on/8 seconds off (begins in the off mode)
- **Flow rates**
 - (Intermittent) 0 to 13 L/min, preset at 8 L/min (ASTM)
 - (Continuous) 0 to 80 L/min (not adjustable)
- **Large, Easy to Read Gauge**
 - Dual scale (kPa and mmHg)
 - 6.5 cm (2.5 in) gauge face
 - Concave lens cover reduces glare
 - Accuracy ± 5% full scale
- **Modular Components**
 - Simplifies service procedures
 - Low cost replacement parts and kits
- **Durable Design and Construction**
 - Few mechanical parts
 - Rugged, shatter-resistant case



Adult Intermittent Suction Unit
 •0-200 mmHg
 •0-26 kPa

- **CE marked**
- **Ordering Options**

Available in various configurations to include probes for connection to source vacuum. Also available with options for collection devices and overflow safety protection.

Specifications*

Flow Rate

Intermittent: 0-13 L/min (preset at 8 L/min)
Continuous: 0-80 L/min (not adjustable) without fittings at full increase

Timing

3 to 30 seconds
Independent on and off
Preset 15 seconds on/8 seconds off
Begins in the "OFF" mode

Vacuum Range

0 to full range of available vacuum

Gauge Range

0-26.7 kPa (0-200 mmHg) and full vacuum

Gauge Accuracy

± 5% Full Scale Deflection

Dimensions

(H x W x D): 16.8 x 8.9 x 12.2 cm(6.6 x 3.5 x 4.8 in)

Weight

0.781 Kg (1 lb 11oz)

Warranty

Three Year

Applied Standards

ISO 10079-3

* Note: Specifications are nominal, subject to change without notice.

Shipping Information

Shipping weight:

1.1 Kg (2 lbs 7 oz) less fittings and adapters

Package Size:

(H x W x D)
20.3 x 8.9 x 12 cm
(8 x 3.5 x 4.7 in)


Packaging:

One per carton

For more information, contact your Ohio Medical sales specialists or authorized agent.

PARTS CONFIGURATION

67XX - 1253 - 9XX

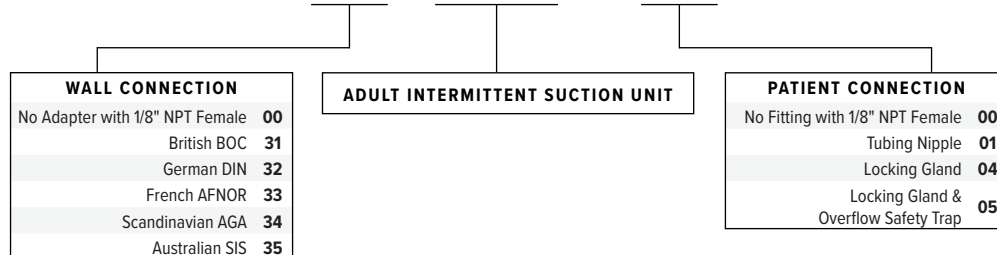

© 2018 Ohio Medical LLC.

This document contains information that is proprietary and confidential to Ohio Medical LLC. Use of this information is under license from Ohio Medical LLC. Any use other than that authorized by Ohio Medical LLC is prohibited.

Ohio Medical and the Ohio Medical logo are registered trademarks of Ohio Medical LLC.



EMERGO EUROPE
Prinsessegracht 20
2514 AP The Hague
The Netherlands

Ohio Medical LLC
1111 Lakeside Drive
Gurnee, IL 60031 USA
866.549.6446 | ohiomedical.com
550041-EN-ISU (Rev. 3) 11/18